
# Adapter Plate 2 (WC-AP2) Multi-Use Plate

Scale 1:1

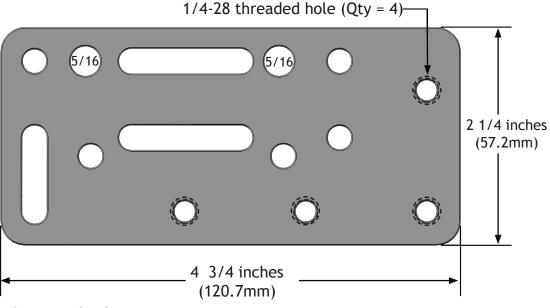

# Set Includes:

(1) Adapter Plate 2 (WC-AP2)

# Hardware - Corpus:

- (4) M5-.8 x 12 Socket Head Cap Screws
- (3) M6-1 x 45 Socket Head Cap Screws
- (2) 0.25 x 0.715 x 0.125 Plastic Washer

## Hardware - Quantum:

(1) 5/16-18 x 1 Hex Cap Screw

#### Hardware - Unitrack:

- (4) M6-1 x 14 Socket Head Cap Screws
- (4) M6-1 x 10 Square Nut


## **Plate to Plate Attachment Options**

- Non-threaded hole to threaded hole = Good
- Non-threaded hole to non-threaded hole with nut = Good
- Threaded hole to threaded hole = Not good

(will damage threads!)

#### **Hole Pairing Options**

- Two holes = **BEST**
- Hole and slot (horizontal or vertical) = Good
- Slot and slot (both horizontal) = Good
- Two vertical slots = NOT GOOD! (could slip!)

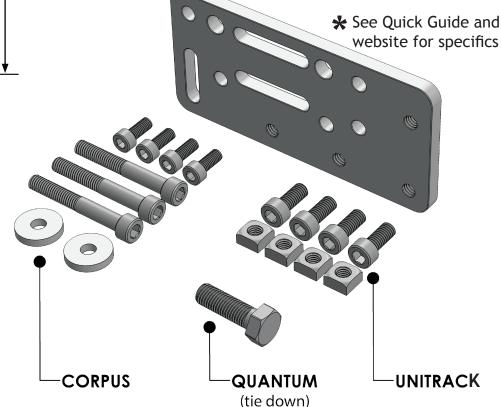

# Mount'n Mover™

by BlueSky Designs

**Weight:** 10.0 oz. (0.28 kg)

**Product:** WC-AP2, Adapter Plate 2

#### **Wheelchair Attachment Overview** \* All plates are considered 'multi-use plates' and can be configured, combined, and used in a variety of ways. **OPTIONAL OFFSET** Wheelchair **Anchor to Chair Bracket** Round tubing \* What you use will Slide track depend on what kind of Hole-spacing chair frame you have. • Tie-Down Adapter Plate 2 T-Nut (WC-AP2) (WC-TN) Or Adapter Plate 4 Adapter Plate 2 (WC-AP4) (WC-AP2) WART . Start with a power or manual Adapter Plate 4 (WC-AP4) Angle Adjust Plate wheelchair (WC-APAP) Round Clamp (WC-R1 or R7/8) Span the bridge L-Angle Extension Plate clamps with a plate (WC-LAE) \* Offset example Adapter Plate 2 (WC-AP2) Bridge Clamps (WC-AP) Plate to span across the (WC-AP4) clamps