


Attaches to:


ECO2, ECOpoint (with eyegaze),
ECO-14, Vanguard, and Pathfinder 2

Device to
device plate

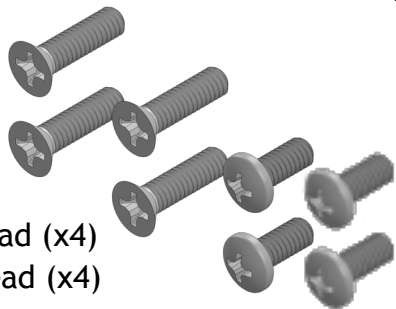


8-32 x 5/8" flathead (x4)
8-32 x 3/8" panhead (x4)

Quick Release
Plate (QRP) to
device plate

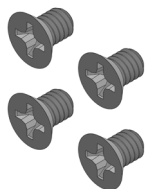

8-32 x 1/4" (x4)

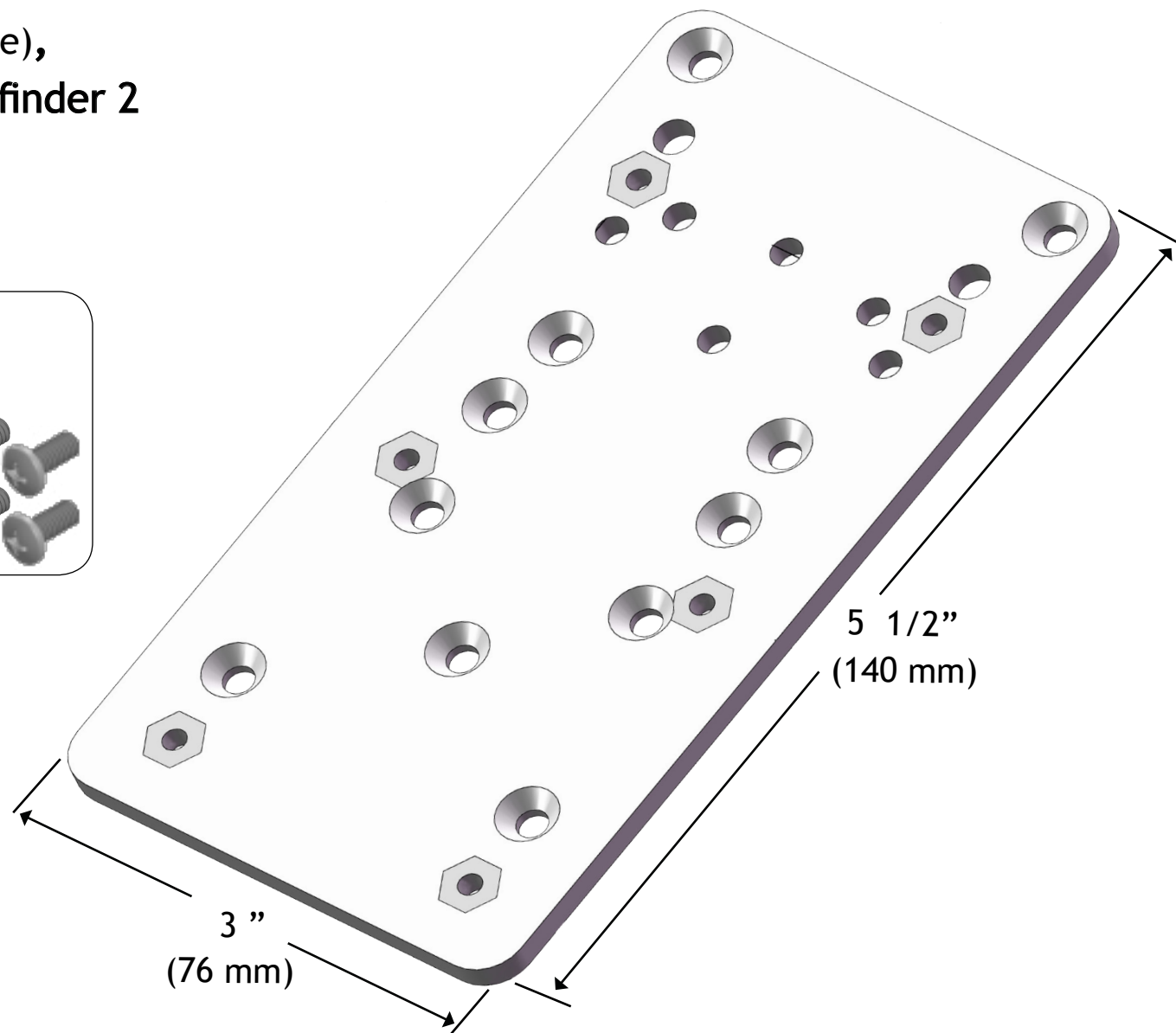

Mount'n Mover™
by BlueSky Designs

mountnmover.com


Product: DP-PRC1, Prentke Romich Co.
Weight: 3.8 oz. (0.11 kg)

1-888-724-7002 (toll free) or 1-612-724-7002

pg 1 of 2

ECO 2 and ECOpoint (with eyegaze)


*For all models, the original mounting plate screws should be removed first

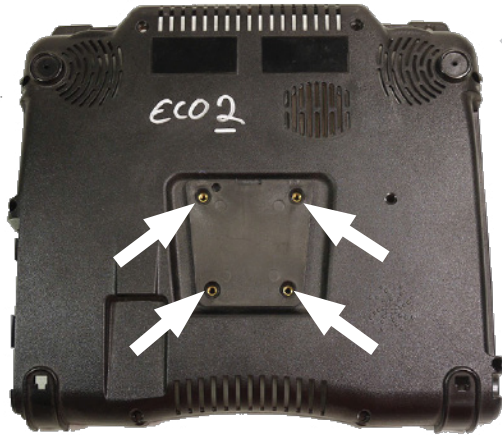


Leave mounting plate in place



Attach device plate

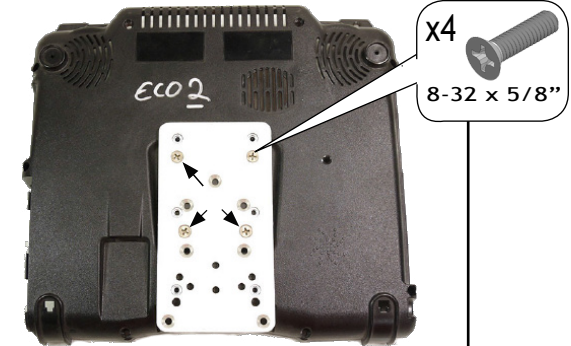

Mount QRP lower when attaching to the **ECOpoint** (eye gaze)

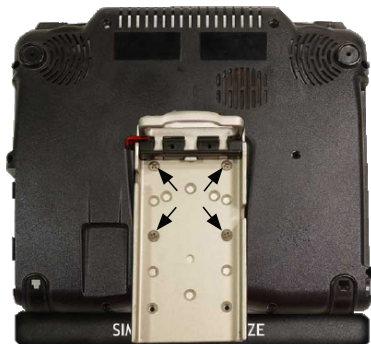

ECOpoint

(Attach before QRP)

Stand+90 or Universal Rotator

ECO2

Attach QRP to ECO2

