



Smartbox GridPad 12 Attachment Instructions

Locking Mount'n Mover (M2-TQ-H) Attachment Instructions

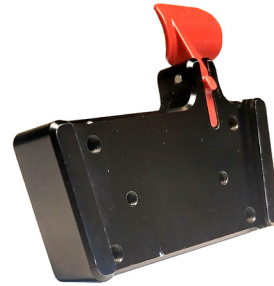

Quick Connect Receiver (QCR)

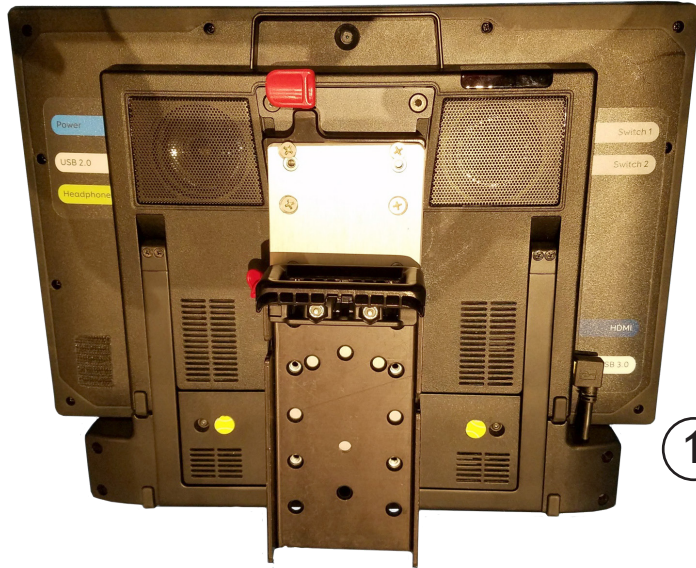

Extension Plate (DP-EXT)




Quick Release Plate (QRP)

- 1 Attach the QCR to the mounting bracket on the back of the GridPad 12
- 2 Use the four screws (8-32 x 1/2") to attach DP-EXT to the QCR (screws included with the QCR).
- 3 Attach the QRP to the threaded inserts on the DP-EXT with four screws (8-32 x 1/4" included with the EXT)

Easy Mover With Quick Connect Attachment Instructions



Easy Mover Dual Arm

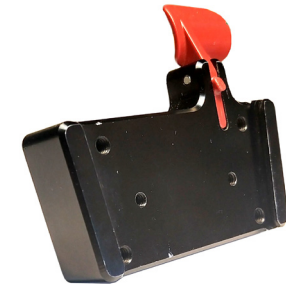

Quick Connect Receiver (QCR)

The Easy Mover with a Quick Connect Receiver attaches directly to the mounting bracket on the back of the GridPad 12.